Multi-lingual Instructions for Use

topdental


Universal Handpiece Oil

Date of revision: 16/01/25

Available in 500ml

UK

Universal Handpiece Oil

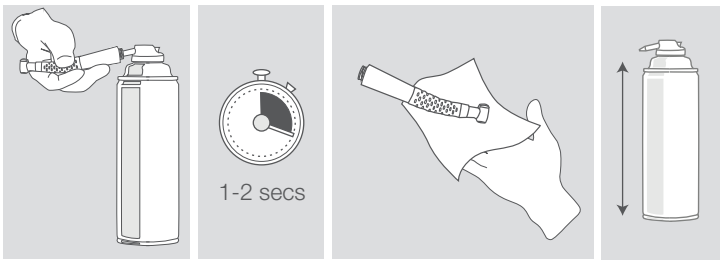All purpose Handpiece Lubricant Spray for maintenance of all dental handpieces including highspeed, slow speed and air motors.

Contains white mineral oil with a trace of peppermint oil.

To use: Insert into the back of the handpiece and spray through for one to two seconds immediately after each use and before sterilisation. Wipe excess thoroughly. Important: Hold can in a vertical position when in use. Refer to handpiece maintenance instructions where applicable. Use only as directed. For professional use only.


Danger. H222: Extremely flammable aerosol. H229: Pressurised container: may burst if heated. P211: Do not spray on an open flame or other ignition source.

Storage. P403+P235: Store in a well ventilated place. Keep cool. P410+P412: Protect from sunlight. Do not expose to temperatures exceeding 50°C/122°F. P102: Keep out of reach of children.

Contraindications. P210: Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking.

Disposal. P251: Do not pierce or burn, even after use.

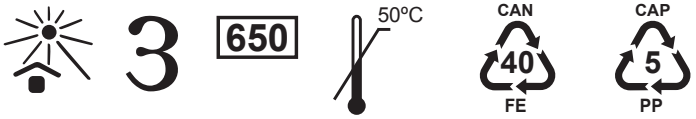UK	**Danger.**
BG	**Опасно.**
CS	**Nebezpečí.**
DA	**Fare.**
DE	**Gefahr.**
EL	**Κίνδυνος.**
ES	**Peligro.**
ET	**Ettevaatust.**
FI	**Vaara.**
FR	**Danger.**
HR	**Opasnost.**
HU	**Veszély.**
IT	**Pericolo.**
LT	**Pavojus.**
LV	**Bīstami.**
NL	**Gevaar.**
NO	**Fare.**
PL	**Niebezpieczeństwo.**
PT	**Perigo.**
RO	**Pericol.**
SK	**Nebezpečenstvo.**
SL	**Nevarno.**
SV	**Fara.**

CE UKCA

Manufactured in the UK by:
Topdental (Products) Ltd.
Unit 1, Holmfield Industrial Estate,
Halifax, W. Yorks, HX2 9TN. UK
Tel +44 (0)1535 652750
www.topdental.co.uk

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