Fig. 1
Abb. 1

Fig. 2
Abb. 2

Fig. 3
Abb. 3

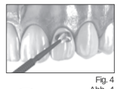
Fig. 4
Abb. 4

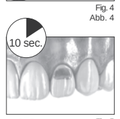
Fig. 5
Abb. 5

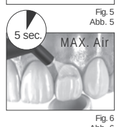
Fig. 6
Abb. 6

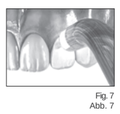
Fig. 7
Abb. 7

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NOT TO BE USED
1. In combination with chemically cured composite resin.
2. In combination with eugenol containing materials as eugenol may hinder G-ænial Bond from setting or bonding properly.
3. In combination with desensitizers as desensitizers may hinder G-ænial Bond from setting or bonding properly.
4. In combination with dual cured luting and core build up composites in case these materials will not be light cured.

DIRECTIONS FOR USE
1. CAVITY PREPARATION
 Prepare tooth using standard techniques. Make sure that the temporary material is thoroughly removed from the surface. Use rubber dam for protection. Dry the prepared tooth surfaces by gently blowing with an air syringe.
 Note:
 For pulp capping, use calcium hydroxide.
2. TECHNIQUE SELECTION
 Select out of the following 2 technique options
 a) Self-etch technique – apply G-ænial Bond to bur cut enamel and dentin, without a separate etching step.
 b) Selective etching of enamel – prior to applying G-ænial Bond to enamel and dentin, etch (bur cut) enamel with 35-40% phosphoric acid gel for 10 seconds, rinse for 5 seconds and gently dry.
 Note:
 Non bur cut enamel should always be treated with a 35-40% phosphoric acid gel for 10 seconds, rinsed with water for 5 seconds and gently dried.
3. APPLICATION
 a) Prior to dispensing, shake the bottle of G-ænial Bond thoroughly (Fig.1). Dispense a few drops into a clean dispensing dish (Fig. 2). Replace bottle cap immediately after use (Fig. 3).
 b) IMMEDIATELY apply to the prepared enamel and dentin surfaces using the disposable applicator (Fig. 4).
 c) Leave undisturbed for 10 seconds after the end of application (Fig.5).
 d) Then, dry thoroughly for 5 seconds with oil free air under MAXIMUM air pressure. Use vacuum suction to prevent splatter of the adhesive (Fig. 6). The final result should be a thin adhesive film with the appearance of frosted glass and which does not visibly move under further air pressure.
 Note:
 1) When removing from refrigerator after long-term storage, leave G-ænial Bond to stand at room temperature for several minutes prior to use.
 2) Apply G-ænial Bond immediately as the material contains a volatile solvent.
 3) Remove excess material remaining on the tooth (other than the surfaces to be bonded) using a sponge or cotton pellet as residue material is hard to remove after light curing.
 4) Should the applied material be contaminated with water, blood or saliva prior to light curing, wash and dry the tooth and repeat the procedure by re-applying material.
4. LIGHT CURING
 Light cure using a visible light curing unit (Fig. 7).
 Irradiation time

Halogen / LED (700 mW/cm²) : 10 seconds
Plasma Arc (2000 mW/cm²) : 3 seconds
G-Light (1200 mW/cm²) : 5 seconds

 In case where the light guide tip is more than 10mm away from the surface to be irradiated, light cure for the following time.

Halogen / LED (700 mW/cm²) : 20 seconds
Plasma Arc (2000 mW/cm²) : 6 seconds
G-Light (1200 mW/cm²) : 10 seconds

 Note:
 1) Light cure completely for an effective bonding strength. Lower light intensity may cause insufficient adhesion.
 2) Use a protective light shield or similar protective eye wear during light curing.

5A. PLACEMENT OF LIGHT CURED COMPOSITES AND COMPOMERS
 After having light cured the adhesive, follow manufacturer's instructions for use in order to place, contour and light cure the composite resin material of choice.

5B. PLACEMENT OF DUAL CURED COMPOSITES
 After having light cured the adhesive, make sure to light cure the dual cured composite material separately. Self curing only will result in insufficient bonding strength.

6. FINISHING
 Adjust restoration, finish and polish using standard techniques.

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CAUTION
1. Wear protective gloves and eye protection.
2. Avoid inhalation or digestion of material.
3. G-ænial Bond is highly flammable liquid and vapour. Do not use near naked flame. Keep away from sources of ignition.
4. G-ænial Bond is volatile. Use in a well ventilated place. Advise patients to breathe through nose.
5. In case of contact with eyes, flush immediately with water and seek medical attention.
6. In case of contact with oral tissue or skin, immediately remove with a sponge or cotton pellet. After the restorative treatment is finished, rinse thoroughly with water.
7. If the tissue contacted by the material turns white or forms a blister, Advise the patient to leave the affected area undisturbed, until the mark disappears, usually in 1-2 weeks. To avoid contact, it is recommended to apply cocoa butter to the area where rubber dam cannot cover.
8. When spilt on the table or floor, immediately wipe off with a dry cloth.
9. Do not mix with other products.
10. Store locked up. Do not store large quantities in one area. Keep away from direct sunlight.
11. Dispose of all waste according to local regulations.

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Prior to use, carefully read the instructions for use. EN

G-ænial® Bond

ONE COMPONENT SELF-ETCHING LIGHT-CURED ADHESIVE

For use only by a dental professional in the recommended indications.

RECOMMENDED INDICATIONS
1. Bonding of light cured composites and acid modified composites (compomers) to tooth structure.
2. Bonding of dual cured luting and core build up composites to tooth structure as long as these materials are light cured.

CONTRAINDICATIONS
1. Pulp capping.
2. In rare cases the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.

NOT TO BE USED
1. In combination with chemically cured composite resin.
2. In combination with eugenol containing materials as eugenol may hinder G-enial Bond from setting or bonding properly.
3. In combination with desensitizers as desensitizers may hinder G-enial Bond from setting or bonding properly.
4. In combination with dual cured luting and core build up composites in case these materials will not be light cured.

DIRECTIONS FOR USE
1. CAVITY PREPARATION
Prepare tooth using standard techniques. Make sure that the temporary material is thoroughly removed from the surface. Use rubber dam for protection. Dry the prepared tooth surfaces by gently blowing with an air syringe.
Note:
For pulp capping, use calcium hydroxide.
2. TECHNIQUE SELECTION
Select out of the following 2 technique options
a) Self-etch technique – apply G-enial Bond to bur cut enamel and dentin, without a separate etching step.
b) Selective etching of enamel – prior to applying G-enial Bond to enamel and dentin, etch (bur cut) enamel with 35-40% phosphoric acid gel for 10 seconds, rinse for 5 seconds and gently dry.
Note:
Non bur cut enamel should always be treated with a 35-40% phosphoric acid gel for 10 seconds, rinsed with water for 5 seconds and gently dried.
3. APPLICATION
a) Prior to dispensing, shake the bottle of G-enial Bond thoroughly (Fig.1). Dispense a few drops into a clean dispensing dish (Fig. 2). Replace bottle cap immediately after use (Fig. 3).
b) IMMEDIATELY apply to the prepared enamel and dentin surfaces using the disposable applicator (Fig. 4).
c) Leave undisturbed for 10 seconds after the end of application (Fig.5).
d) Then, dry thoroughly for 5 seconds with oil free air under MAXIMUM air pressure. Use vacuum suction to prevent splatter of the adhesive (Fig. 6). The final result should be a thin adhesive film with the appearance of frosted glass and which does not visibly move under the air pressure.
Note:
1) When removing from refrigerator after long-term storage, leave G-enial Bond to stand at room temperature for several minutes prior to use.
2) Apply G-enial Bond immediately as the material contains a volatile solvent.
3) Remove excess material remaining on the tooth (other than the surfaces to be bonded) using a sponge or cotton pellet as residue material is hard to remove after light curing.
4) Should the applicator material be contaminated with water, blood or saliva prior to light curing, wash and dry the tooth and repeat the procedure by re-applying material.
4. LIGHT CURING
Light cure using a visible light curing unit (Fig. 7).
Irradiation time

Halogen / LED (700 mW/cm²) : 10 seconds
Plasma Arc (2000 mW/cm²) : 3 seconds
G-Light (1200 mW/cm²) : 5 seconds

In case where the light guide tip is more than 10mm away from the surface to be irradiated, light cure for the following time.

Halogen / LED (700 mW/cm²) : 20 seconds
Plasma Arc (2000 mW/cm²) : 6 seconds
G-Light (1200 mW/cm²) : 10 seconds

Note:
1) Light cure completely for an effective bonding strength. Lower light intensity may cause insufficient adhesion.
2) Use a protective light shield or similar protective eye wear during light curing.
5A. PLACEMENT OF LIGHT CURED COMPOSITES AND COMPOMERS
After having light cured the adhesive, follow manufacturer's instructions for use in order to place, contour and light cure the composite resin material of choice.
5B. PLACEMENT OF DUAL CURED COMPOSITES
After having light cured the adhesive, make sure to light cure the dual cured composite material separately. Self curing only will result in insufficient bonding strength.
6. FINISHING
Adjust restoration, finish and polish using standard techniques.

STORAGE
Store at room temperature (1-28°C) (33.8-82.4°F).
If not in use for a prolonged period of time, store in refrigerator.
Shelf life: 2 years from date of manufacture.

PACKAGE
1. G-enial Bond Bottle Kit:
5 mL liquid (1), disposable dispensing dish (20), disposable applicator fine (50)
2. G-enial Bond Bottle Refill:
5 mL liquid (1)
3. G-enial Bond 3-bottle Pack:
5 mL liquid (3)

CAUTION
1. Wear protective gloves and eye protection.
2. Avoid inhalation or digestion of material.
3. G-enial Bond is highly flammable liquid and vapour. Do not use near naked flame. Keep away from sources of ignition.
4. G-enial Bond is volatile. Use in a well ventilated place. Advise patients to breathe through nose.
5. In case of contact with eyes, flush immediately with water and seek medical attention.
6. In case of contact with oral tissue or skin, immediately remove with a sponge or cotton pellet. After the restorative treatment is finished, rinse thoroughly with water.
7. If the tissue contacted by the material turns white or forms a blister, advise the patient to leave the affected area undisturbed, until the mark disappears, usually in 1-2 weeks. To avoid contact, it is recommended to apply cocoa butter to the area where rubber dam cannot cover.
8. When spilt on the table or floor, immediately wipe off with a dry cloth.
9. Do not mix with other products.
10. Store locked up. Do not store large quantities in one area. Keep away from direct sunlight.
11. Dispose of all waste according to local regulations.

Some products referenced in the present IFU may be classified as hazardous according to GHS.
Always familiarize yourself with the safety data sheets available at:
http://www.gceurope.com/download
They can also be obtained from your supplier.

Last revised : 08/2015

CE 0086

MANUFACTURED by
GC CORPORATION
76-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8585, Japan

DISTRIBUTED by
GC CORPORATION
76-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8585, Japan

EU: GC EUROPE N.V.
Researchpark Haasrode-Leuven 1240, Interleuvenlaan 33,
B-3001 Leuven, Belgium TEL: +32 16 74 10 00

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