COMPONENTS
КОМПОНЕНТИ
KOMPONENTY
DIJELOVI
RÉSZEK
CZĘŚCI
COMPONENTE
КОМПОНЕНТЫ

① Capsule Cover (polyacetal)
Капак за капсулата (полиацетал)
Povrch kapsle (polyacetál)
Poklopac kapsule (poliacetal)
Kapszula fedő (poliacetál)
Osłona kapsułki (poliacetal)
Capac capsulă (poliacetal)
Заглушка капсулы (полиоксиметилен)
② Plunger (stainless steel)
Бутало (стомана)
Píst (nerezová ocel)
Klip (nehrdjúci čelik)
Dugattyú (rozdamentesacél)
Tlok (stal nierdzewna)
Piston (oțel inoxidabil)
Поршень (нержавеющая сталь)
③ Capsule holder (stainless steel)
Държач на капсулите (стомана)
Držák kapslí (nerezováocel)
Držiak kapsule (nehrdjúci čelik)
Kapszula tartó (rozdamentes acél)
Uchwyt na kapsułki (stal nierdzewna)
Suport capsulă (oțel inoxidabil)
Фиксатор капсулы (нержавеющая сталь)
④ Release button (polyacetal)
Освобождаващ бутон (полиацетал)
Uvolňovací tlačítko (polyacetál)
Poluga za otpuštanje (poliacetal)
Kioldó gomb (poliacetal)
Przycisk zwalniający (poliacetal)
Buton declanșator (poliacetal)
Освобождающая кнопка (полиоксиметилен)
⑤ Grip (stainless steel)
Дръжка (стомана)
Rukojeť (nerezová ocel)
Ručka (nehrdjúci čelik)
Tartó kar (rozdamentes acél)
Uchwyt (stal nierdzewna)
Mâner (oțel inoxidabil)
Рукоять (нержавеющая сталь)
⑥ Lever (stainless steel)
Лост (стомана)
Páka (nerezová ocel)
Poluga (nehrdjúci čelik)
Adagoló kar (rozdamentes acél)
Dźwignia (stal nierdzewna)
Manetă (oțel inoxidabil)
Рычаг (нержавеющая сталь)

Fig. 1 / Фиг. 1/ obr. 1 / slika 1/ 1. ábra / ryc. 1 / Рис. 1

Fig. 2 / Фиг. 2 / obr. 2 / slika 2 / 2. ábra / ryc. 2 / Рис. 2

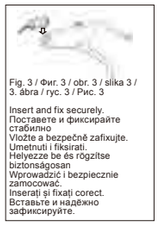
Fig. 3 / Фиг. 3 / obr. 3 / slika 3 / 3. ábra / ryc. 3 / Рис. 3

Insert and fix securely.
Поставете и фиксирайте стабилно
Vložte a bezpečně zafixujte.
Umetnuti i fiksirati.
Helyezze be és rögzítse biztonságosan
Wprowadzić i bezpiecznie zamocować.
Inserați și fixați corect.
Вставьте и надёжно зафиксируйте.

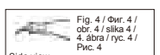
Fig. 4 / Фиг. 4 / obr. 4 / slika 4 / 4. ábra / ryc. 4 / Рис. 4

Side view
Страничен изглед
Boční pohled
Izgled sa strane
Oldalnézet
Widok z boku
Vedere laterală
Вид сбоку

Fig. 5 / Фиг. 5 / obr. 5 / slika 5 / 5. ábra / ryc. 5 / Рис. 5

Top view
Изглед отгоре
Pohled shora
Izgled odozgo
Felülnézet
Widok z góry
Vedere în plan
Вид сверху

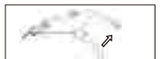
Fig. 6 / Фиг. 6 / obr. 6 / slika 6 / 6. ábra / ryc. 6 / Рис. 6

2 priming strokes + 2 slow extruding strokes.
2 първични щраквания + 2 бавни екструдиращи щраквания.
2 vstřikovací tahy + 2 pomalé vytlačovací tahy.
2 pritiska za aktiviranje + 2 polagana pritiska za istiskivanje
Két beélesítő nyomás majd két lassú nyomás a kiadagoláshoz.
2 pociągnięcia do aktywowania + 2 pociągnięcia do powolnego wytłaczania.
2 apăsări pentru amorsare + 2 apăsări ușoare pentru dozare.
2 нажатия для фиксации + 2 медленных нажатия для извлечения.

Fig. 7 / Фиг. 7 / obr. 7 / slika 7 / 7. ábra / ryc. 7 / Рис. 7

① Retract the plunger.
Връщане на буталото
Zatáhněte píst.
Uvučí klip.
Húzza vissza a dugattyút.
Cofnąć tłok.
Retragerea pistonului.
Отведите назад поршень
② Pull the capsule.
Издърпайте капсулата
Vytáhněte kapsli.
Povučí kapsulu.
Wyciągnąć kapsułkę
Scoaterea capsulei
Извлеките капсулу

GC CAPSULE APPLIER Ⅳ

Prior to use, carefully read the instructions for use. EN

For use only by a dental professional in the indications for use.

USE
Dispensing of encapsulated GC cements.

DIRECTIONS FOR USE
1. Capsule activation
 Before activation, shake the capsule or tap its side on a hard surface to loosen the powder (Fig. 1).
 To activate the capsule, push the plunger until it is fully depressed with the main body and hold it down for 2 seconds (Fig. 2).
2. Dispensing of mixed cement
 a) After mixing in a mixer or amalgamator, insert and fix the capsule securely in the capsule holder (Fig. 3, 4, 5).
 NOTE:
 The capsule cannot be inserted if the plunger is extended. Retract the plunger by pressing the release button.
 b) Make 2 strokes with lever to prime the capsule.
 Make the third stroke slowly to check if the cement is extruded in the nozzle. Direct the capsule nozzle into the restoration or prepared tooth cavity and make slowly the fourth stroke to allow smooth dispensing of the cement (Fig. 6).
 Avoid putting too much force onto the applier.
 NOTE:
 1) To adjust the direction of the nozzle, hold the applier with the capsule towards you and turn the capsule body.
 2) In case of difficulty in extrusion, remove the capsule and check the applier. If the applier is in good condition, replace with a new capsule.
3. Removal of the capsule
 Press the release button to retract the plunger. Remove the used capsule from the applier by twisting and pulling upwards (Fig. 7).

CLEANING / DISINFECTING
Remove any remaining cement immediately from the applier. Clean the applier with a gauze soaked in alcohol or in an ultrasonic bath. In case of using chemical disinfectant / sterilizing agent, follow the manufacturer's instructions carefully.
NOTE:
When immersing the applier in the chemical disinfectant for a prolonged period of time, it is recommended to use a disinfectant solution which contains rust inhibitor.

STERILIZATION
The applier can be autoclaved at 121-134°C (250-273°F) for 10 to 30 minutes. Please follow manufacturer's instructions.
NOTE:
The operating (sterilizing and drying) temperature should not exceed 134°C (273°F), as this may cause deformation of the polyacetal parts.

CAUTION
1. Prior to use, sterilize the applier.
2. Take care not to get your hand or fingers caught between the lever and the grip.
3. In rare cases, the product may cause sensitivity in some people. If any such reactions are experienced, discontinue the use of the product and refer to a physician.
4. Personal protective equipment (PPE) such as gloves, face masks and safety eyewear should always be worn.

STORAGE
Keep GC CAPSULE APPLIER Ⅳ in a dry place away from high moisture to prevent from rusting.

PACKAGE
GC CAPSULE APPLIER Ⅳ - 1 piece

Undesired effects- Reporting:
If you become aware of any kind of undesired effect, reaction or similar events experienced by use of this product, including those not listed in this instruction for use, please report them directly through the relevant vigilance system, by selecting the proper authority of your country accessible through the following link:
https://ec.europa.eu/growth/sectors/medical-devices/contacts_en
as well as to our internal vigilance system:
vigilance@gc.dental
In this way you will contribute to improve the safety of this product.

Last revised: 07/2024

DISTRIBUTED by
GC CORPORATION
76-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8585, Japan

EU : GC EUROPE N.V.
Researchpark Haasrode-Leuven 1240, Interleuvenlaan 33, B-3001 Leuven, Belgium TEL: +32 16 74 10 00

MANUFACTURED by
MCC Advanced Moldings Co., Ltd.
1-1, Marunouchi 1-chome, Chiyoda-ku, Tokyo, Japan

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