Multi-lingual Instructions for Use

topdental

Alcohol Free Surface Disinfectant

Date of revision: 10/09/24

Available in 5 Litre, 2.5L and 1 Litre / 500ml trigger spray bottles

UK

Alcohol-Free Surface Disinfectant

MD **For the surface disinfection of medical devices, including but not limited to, dental chairs, examination lights and patient transport.** Not intended for use on medical and dental instruments. Specially formulated for use in the medical and dental environment.

BPR **For the disinfection of surfaces.** Use biocides safely and sustainably. This product should be used in accordance with the product label. PT2 Formulation type: AL.

To use: Spray clean surfaces and objects to be disinfected from a distance of 30cm so that they are completely moistened and let the disinfectant take action over 60 seconds (to achieve full efficacy leave for 5 minutes). Allow surface moisture to dry. Wipe with clean tissue if necessary. Contains detergent. For professional use only.


Efficacy: Bactericidal (5mins). Yeasticidal (5 mins). Virucidal activity against enveloped viruses (60 secs): Vaccinia virus, Hepatitis B virus (HBV), Hepatitis C virus (HCV), Hepatitis Delta virus, HIV, Influenza virus, Coronavirus, Herpesviridae, Filoviridae (Ebola), Rubella, Rabies, Measles, Poxviridae.


Tested to standards: EN16777, EN14885, EN14476, EN13727, EN13697, EN13624.


Composition: Contains Aqua, Alkyl (C12-16) dimethylbenzyl ammonium chloride (ADBAC/BKC (C12-16)) (68424-85-1) 0.48% w/w (4.8g/kg), Didecyldimethylammonium Chloride (7173-51-5) 0.22% w/w (2.2g/kg), Perfumes (contain Amyl cinnamal, Limonene, Alpha-Isomethyl ionone).

Handling and Storage. Ensure there is sufficient ventilation of the area. The wearing of protective gloves is recommended. P403+P233: Store in a well-ventilated place. Keep container tightly closed. P235: Keep cool. **Contraindications.** Mild steel, aluminium, copper. **Disposal.** Dispose of contents/container in accordance with local regulations. Do not discharge liquid into drains or rivers. UK number for poison advice: 111. Poisons Information Ireland: For information or to report a poisoning incident contact The National Poisons Information Centre (01 8092166).

H412: Harmful to aquatic life with long lasting effects. P273: Avoid release to the environment. P301+P310: IF SWALLOWED: Immediately call a doctor. P305+P351+P338: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

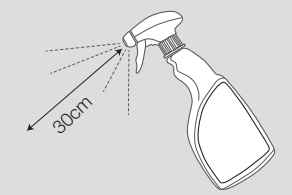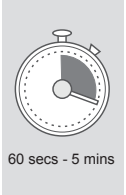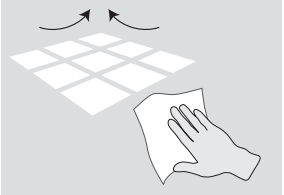Efficacy

Microorganism	Time
Bactericidal	**5 min**
Pseudomonas aeruginosa	5 min
Staphyloccocus aureus	5 min
Enterococcus hirae	5 min
Yeasticidal	**5 min**
Candida albicans	5 min
Enveloped viruses	**60 sec**
Vaccinia virus	60 sec
Hepatitis B virus (HBV)	60 sec
Hepatitis C virus (HCV)	60 sec
Hepatitis Delta virus	60 sec
HIV	60 sec
Influenza virus	60 sec
Coronavirus	60 sec
Herpesviridae	60 sec
Filoviridae (Ebola)	60 sec
Rubella virus	60 sec
Rabies virus	60 sec
Measles virus	60 sec
Poxviridae	60 sec

5L Refill size may only be used to refill dispensers which are labelled with a DAFM approved label for the product.

CE 1639

UK CA 0120

Manufactured in the UK by:
Topdental (Products) Ltd.
Unit 1, Holmfield Industrial Estate,
Halifax, W. Yorks, HX2 9TN. UK
Tel +44 (0)1535 652750
www.topdental.co.uk

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