Multi-lingual Instructions for Use

topdental

Alcohol Free Surface Disinfectant Wipes

Date of revision: 12/09/24

Available in 200 and 160 wipe dispenser tubs and refill bags

UK


Alcohol-Free Surface Disinfectant Wipes

MD

For the surface disinfection of medical devices, including but not limited to, dental chairs, examination lights and patient transport. Not intended for use on medical and dental instruments. Specially formulated for use in the medical and dental environment.

BPR

For the disinfection of surfaces. Use biocides safely and sustainably. This product should be used in accordance with the product label. PT2 Formulation type: AL.

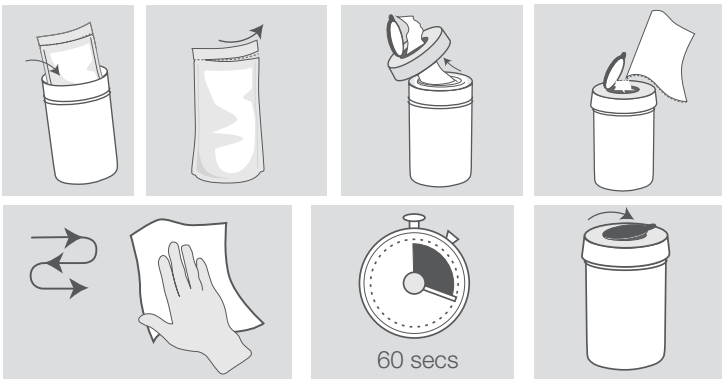To use: 1. Place refill bag in dispenser tub, 2. Cut open air tight plastic pouch across the top of the pouch, 3. Feed the centre wipe through the hole in the cap and replace, 4. Pull wipe at a 45º angle. Using sufficient wipes for each application wipe the clean surfaces to be disinfected so that they are completely moistened and let the disinfectant take action over 60 seconds. Allow surface moisture to dry or wipe with clean tissue. Contains detergent. For professional use only.

Efficacy (60 secs): Bactericidal. Yeasticidal. Virucidal activity against enveloped viruses: Vaccinia virus, Hepatitis B virus (HBV), Hepatitis C virus (HCV), Hepatitis Delta virus (HDV), Human Immunodeficiency Virus (HIV), Influenza virus, Coronavirus, Herpesviridae, Filoviridae (Ebola), Rubella, Rabies, Measles, Poxviridae.

Tested to Standards: EN16777, EN16615, EN14885, EN14476, EN13727, EN13624.

Composition: Contains Aqua, Alkyl (C12-16) dimethylbenzyl ammonium chloride (ADBAC/BKC (C12-16)) (68424-85-1) 0.48% w/w (4.8g/kg), Didecyldimethylammonium Chloride (7173-51-5) 0.22% w/w (2.2g/kg), Perfumes (contain Amyl cinnamal, Limonene, Alpha-Isomethyl ionone).

Handling and Storage. Ensure there is sufficient ventilation of the area. The wearing of protective gloves is recommended. P403+P233: Store in a well-ventilated place. Keep container tightly closed. P235: Keep cool. **Contraindications.** Mild steel, aluminium, copper. **Disposal.** Dispose of contents/container in accordance with local regulations. Do not discharge liquid into drains or rivers. UK number for poison advice: 111. Poisons Information Ireland: For information or to report a poisoning incident contact The National Poisons Information Centre (01 8092166). Do not flush. Non maceratable.

H412: Harmful to aquatic life with long lasting effects. P273: Avoid release to the environment. P301+P310: IF SWALLOWED: Immediately call a doctor. P305+P351+P338: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing.

Efficacy

Microorganism	Time
Bactericidal	**60 sec**
Pseudomonas aeruginosa	60 sec
Staphyloccocus aureus	60 sec
Enterococcus hirae	60 sec
Yeasticidal	**60 sec**
Candida albicans	60 sec
Enveloped viruses	**60 sec**
Vaccinia virus	60 sec
Hepatitis B virus (HBV)	60 sec
Hepatitis C virus (HCV)	60 sec
Hepatitis Delta virus	60 sec
HIV	60 sec
Influenza virus	60 sec
Coronavirus	60 sec
Herpesviridae	60 sec
Filoviridae (Ebola)	60 sec
Rubella virus	60 sec
Rabies virus	60 sec
Measles virus	60 sec
Poxviridae	60 sec

CE 1639

UK CA 0120

Manufactured in the UK by:
Topdental (Products) Ltd.
Unit 1, Holmfield Industrial Estate,
Halifax, W. Yorks, HX2 9TN. UK
Tel +44 (0)1535 652750
www.topdental.co.uk

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