Multi-lingual Instructions for Use

topdental


Surface Disinfectant Wipes - Alcohol

Date of revision: 28/08/24

Available in 200 and 160 wipe dispenser tubs and refill bags

UK

Surface Disinfectant Wipes

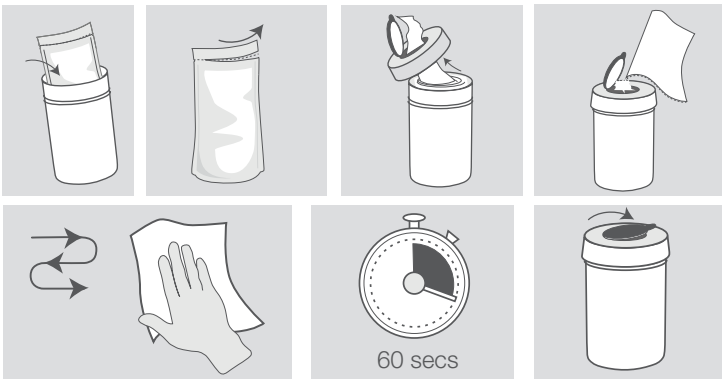MD **For the surface disinfection of non-invasive medical devices.** Alcohol is not recommended for use on dental/medical instruments.

BPR **For the disinfection of surfaces.** Use biocides safely and sustainably. This product should be used in accordance with the product label. PT2 Formulation type: AL.

To use: 1. Place refill bag in dispenser tub, 2. Cut open air tight plastic pouch across the top of the pouch, 3. Feed the centre wipe through the hole in the cap and replace, 4. Pull wipe at a 45º angle. Using sufficient wipes for each application wipe the clean surfaces to be disinfected so that they are completely moistened and let the disinfectant take action for 60 seconds. Allow surface moisture to evaporate. For professional use only.

Efficacy (60 secs): Bactericidal including MRSA, VRE. Yeasticidal. Mycobactericidal. Tuberculocidal. Viruses to include: Non-enveloped viruses: Murine norovirus, Adenovirus 5. Enveloped viruses: Hepatitis B virus (HBV), Hepatitis C virus (HCV), Hepatitis Delta virus, HIV, Influenza virus, Coronavirus, Herpesviridae, Filoviridae (Ebola), Rubella, Rabies, Measles, Poxviridae.

Tested to Standards: EN16615, EN14885, EN14476, EN14348, EN13727, EN13624.

Composition: Ethanol (64-17-5) 63% w/w (630g/kg), Aqua, Perfumes (Contains Linalool, Citral, Geraniol, Alpha-Isomethyl ionone).

Contraindications. Do not use on surfaces sensitive to alcohol. **Disposal.** Dispose of contents/container in accordance with local regulations. Do not discharge into drains or rivers. Do not flush. Non maceratable. UK number for poison advice: 111. Poisons Information Ireland: For information or to report a poisoning incident contact The National Poisons Information Centre (01 8092166).

Danger. H225: Highly flammable liquid and vapour. H319: Causes serious eye irritation. P280: Wear eye protection, protective gloves. P337+P313: If eye irritation persists: Get medical advice/attention. P370+P378: In case of fire: Use alcohol resistant foam to extinguish. **Storage.** P210: Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. P403+P233: Store in a well-ventilated place. Keep container tightly closed. P235: Keep cool.

Efficacy

Microorganism	Time
Bactericidal	**60 sec**
Pseudomonas aeruginosa	60 sec
Staphyloccocus aureus	60 sec
Enterococcus hirae	60 sec
E. coli	60 sec
MRSA	60 sec
VRE	60 sec
Yeasticidal	**60 sec**
Candida albicans	60 sec
Enveloped viruses	**60 sec**
Hepatitis B virus (HBV)	60 sec
Hepatitis C virus (HCV)	60 sec
HIV	60 sec
Measles virus	60 sec
Hepatitis Delta virus	60 sec
Influenza virus	60 sec
Coronavirus	60 sec
Herpesviridae	60 sec
Filoviridae (Ebola)	60 sec
Rubella virus	60 sec
Rabies virus	60 sec
Poxviridae	60 sec
Non enveloped viruses	**60 sec**
Adenovirus 5	60 sec
Murine norovirus	60 sec
Mycobactericidal	**60 sec**
Mycobacterium avium	60 sec
Tuberculocidal	**60 sec**
Mycobacterium terrae	60 sec

UK	Danger.
BG	Опасно.
CS	Nebezpečí.
DA	Fare.
DE	Gefahr.
EL	Κίνδυνος.
ES	Peligro.
ET	Ettevaatust.
FI	Vaara.
FR	Danger.
HR	Opasnost.
HU	Veszély.
IT	Pericolo.
LT	Pavojus.
LV	Bīstami.
NL	Gevaar.
NO	Fare.
PL	Niebezpieczeństwo.
PT	Perigo.
RO	Pericol.
SK	Nebezpečenstvo.
SL	Nevarno.
SV	Fara.

CE 1639

UK CA 0120

Manufactured in the UK by:
Topdental (Products) Ltd.
Unit 1, Holmfield Industrial Estate,
Halifax, W. Yorks, HX2 9TN. UK
Tel +44 (0)1535 652750
www.topdental.co.uk

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