Multi-lingual Instructions for Use

topdental


Instrument Disinfectant and Cleaner


Date of revision: 20/01/25


Available in 5 Litre / 2.5 Litre / 1 Litre dispenser bottle

UK

Instrument Disinfectant and Cleaner

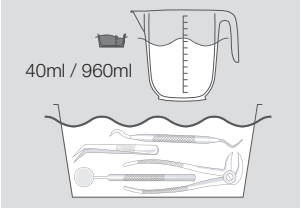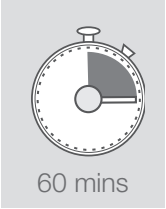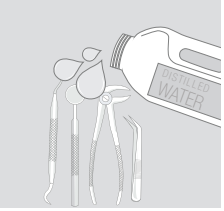Disinfecting and cleaning concentrate for instruments, burs and endodontic instruments and items. Chloride, aldehyde and phenol free with a wide material compatibility and broad microbiological efficacy.

To Use: Add 20ml ± 1ml to 980ml of clean water for a 2% solution, 40ml ± 1ml to 960ml for a 4% solution. **To achieve full efficacy use a 4% solution for a 60 minute soaking time.** After disinfection, rinse instruments thoroughly with RO or distilled water prior to further processing and sterilisation. Follow the reprocessing instructions from the instrument manufacturer. Change the solution daily. For professional use only.

Efficacy 2%: Bactericidal including MRSA (15 min). Yeasticidal (15 min). All enveloped viruses (15 min) to include: Vaccinia virus, HBV, HCV, HIV, Hepatitis Delta virus, Influenza virus, Coronavirus, Herpesviridae, Filoviridae (Ebola), Rubella, Rabies, Measles, Poxviridae.

Efficacy 4%: *As above and to include* Non-enveloped viruses (30 min): Murine norovirus, Adenovirus 5. Mycobactericidal & Tuberculocidal (60 min).

Tested to Standards: EN17111, EN14885, EN14561, EN14562, EN14563, EN14476, EN14348, EN13727, EN13624.

Composition: Aqua, N-(3-Aminopropyl)-N-Dodecylpropane-1,3-Diamine (4.1% w/w 2372-82-9), DMPAP (3.5% w/w), Alcohols, C9-11, Ethoxylated, Tetrasodium N,N-Bis(Carboxylatomethyl)-L-Glutamate, Monoethanolamine, Corrosion Inhibitor, Perfume contains CINEOLE - Cineols - Eucolytus; HYDRATROPALDEHYDE; TRANS-MENTHONE; p-menth-1-en-4-ol., Dye.

Storage. Store in a well-ventilated place. Keep container tightly closed. Keep cool. **Contraindications.** Heat. **Disposal.** Do not discharge concentrate into drains or rivers.

Danger. H315: Causes skin irritation. H318: Causes serious eye damage. H410: Very toxic to aquatic life with long lasting effects. P264: Wash hands thoroughly after handling. P280: Wear protective gloves, eye protection, protective clothing. P305+P351+P338+P310: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a doctor. P332+P313: If skin irritation occurs: Get medical advice/ attention. P362+P364: Take off contaminated clothing and wash it before reuse. P391: Collect spillage.

Efficacy

Microorganism	Time	Conc.
Bactericidal	**15 min**	**2%**
Pseudomonas aeruginosa	15 min	2%
Staphyloccocus aureus	15 min	2%
Enterococcus hirae	15 min	2%
E. coli	15 min	2%
Salmonella typhimurium	15 min	2%
Legionella pneumophila	15 min	2%
MRSA	15 min	2%
Yeasticidal	**15 min**	**2%**
Candida albicans	15 min	2%
Enveloped viruses	**15 min**	**2%**
Vaccinia virus	15 min	2%
Hepatitis B virus (HBV)	15 min	2%
Hepatitis C virus (HCV)	15 min	2%
HIV	15 min	2%
Measles virus	15 min	2%
Hepatitis Delta virus	15 min	2%
Influenza virus	15 min	2%
Coronavirus	15 min	2%
Herpesviridae	15 min	2%
Filoviridae (Ebola)	15 min	2%
Rubella virus	15 min	2%
Rabies virus	15 min	2%
Non enveloped viruses	**30 min**	**4%**
Adenovirus 5	30 min	4%
Murine norovirus	30 min	4%
Mycobactericidal	**60 min**	**4%**
Mycobacterium avium	60 min	4%
Tuberculocidal	**60 min**	**4%**
Mycobacterium terrae	60 min	4%

UK	Danger.
BG	Опасно.
CS	Nebezpečí.
DA	Fare.
DE	Gefahr.
EL	Κίνδυνος.
ES	Peligro.
ET	Ettevaatust.
FI	Vaara.
FR	Danger.
HR	Opasnost.
HU	Veszély.
IT	Pericolo.
LT	Pavojus.
LV	Bīstami.
NL	Gevaar.
NO	Fare.
PL	Niebezpieczeństwo.
PT	Perigo.
RO	Pericol.
SK	Nebezpečenstvo.
SL	Nevarno.
SV	Fara.

CE 1639

UK CA 0120

Manufactured in the UK by:
Topdental (Products) Ltd.
Unit 1, Holmfield Industrial Estate,
Halifax, W. Yorks, HX2 9TN. UK
Tel +44 (0)1535 652750
www.topdental.co.uk

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