Multi-lingual Instructions for Use

topdental


Ultrasonic Cleaner

Date of revision: 09/08/24

Available in 5 Litre / 2.5 Litre / 1 Litre dispenser bottle

UK

Ultrasonic Cleaner

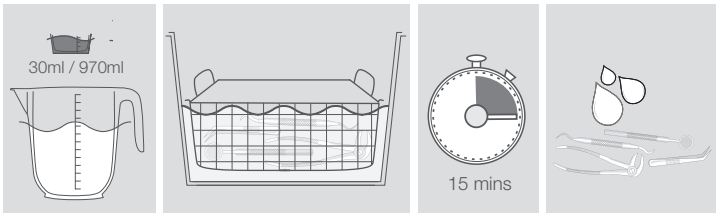Enzyme based concentrate for the pre-cleaning and decontamination of the non-porous surfaces of hand instruments and forceps, rotary instruments, mirrors etc. Designed for use in the Ultrasonic bath. Specially formulated for use in the Medical and Dental environment.

To Use: Add 30ml ± 1ml of the concentrate to 970ml of tap water in the ultrasonic bath. (3% solution). Immerse the instruments ensuring that all surfaces are covered by the solution, add hollow objects in such a way that air escapes. Run the ultrasonic cleaner for a minimum of 15 minutes and according to the manufacturers instructions. After immersion, remove the instruments and rinse them thoroughly prior to terminal sterilisation or high level disinfection. Prepare a fresh solution daily or sooner if solution is very dirty. For professional use only.

Composition: Aqua, Propan-2-ol, Didecyldimethylammonium chloride, Sodium n-lauroylsarcosinate, Subtilisin, Dye.

Storage. P210: Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. P403+P233: Store in a well-ventilated place. Keep container tightly closed. P235: Keep cool. **Contraindications.** Aluminium, copper. Not recommended for plastics and rubbers. **Disposal.** Do not discharge concentrate into drains or rivers.

Warning. H226: Flammable liquid and vapour. H315: Causes skin irritation. H319: Causes serious eye irritation. H336: May cause drowsiness or dizziness. P302+350: IF ON SKIN: Gently wash with plenty of soap and water. P305+351+338: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. P280: Wear protective gloves/protective clothing/eye protection/face protection. P261: Avoid breathing dust/fumes/gas/mist/vapours/spray. P304+340: IF INHALED: Remove person to fresh air and keep comfortable for breathing.

UK	**Warning.**
BG	**Внимание.**
CS	**Varování.**
DA	**Warnung.**
DE	**Achtung.**
EL	**Προσοχή.**
ES	**Atención.**
ET	**Hoiatus.**
FI	**Varoitus.**
FR	**Attention.**
HR	**Upozorenje.**
HU	**Figyelem.**
IT	**Attenzione.**
LT	**Įspėjimas.**
LV	**Brīdinājums.**
NL	**Waarschuwing.**
NO	**Advarsel.**
PL	**Uwaga.**
PT	**Atenção.**
RO	**Atenţie.**
SK	**Pozor.**
SL	**Pozor.**
SV	**Varning.**

CE UKCA

Manufactured in the UK by:
Topdental (Products) Ltd.
Unit 1, Holmfield Industrial Estate,
Halifax, W. Yorks, HX2 9TN. UK
Tel +44 (0)1535 652750
www.topdental.co.uk

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