Multi-lingual Instructions for Use

topdental


Impression Disinfectant Powder


Date of revision: 17/10/24


Available in 700g


UK

Impression Disinfectant Powder

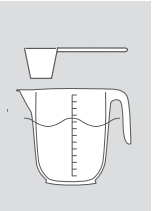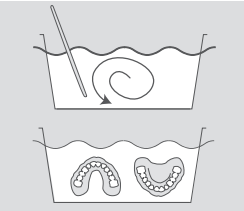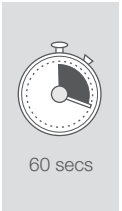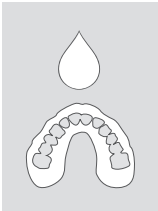Fast acting disinfectant for all kinds of impression material such as silicone, alginates and polyether.

Directions for use: Dissolve 1 level measure of powder (20g ± 2g) in 980ml of warm tap water (30 - 40°C) and mix for 60 seconds. Place the impressions into the solution and fully immerse for the 60 second contact time. Do not exceed this recommended time. Rinse immediately with copious amounts of running water. Change solution daily. For professional use only.

Efficacy (60 secs): Bactericidal including MRSA. Yeasticidal. Viruses to include: Non-enveloped viruses: Murine norovirus, Adenovirus 5. Enveloped viruses: Hepatitis B virus (HBV), Hepatitis C virus (HCV), Hepatitis Delta virus, HIV, Influenza virus, Coronavirus, Herpesviridae, Filoviridae (Ebola), Rubella, Rabies, Measles. Poxviridae.

Tested to Standards: EN 17111, EN14885, EN14476, EN13727, EN13697, EN13624.

Composition: TAED, Sodium Percarbonate, Citric Acid Monohydrate, Sodium Alkyl Aryl Sulphonate, Sodium Sulphate Anhydrous, Tetrasodium (1-hydroxyethylidene) bisphosphonate, Sodium Nitrite.

Contraindications. Aluminium, Copper. **Disposal.** Do not discharge into drains or rivers.

Danger. H318: Causes serious eye damage. H272: May intensify fire; oxidiser. H411: Toxic to aquatic life with long lasting effects. P280: Wear eye protection, protective gloves. P305+P351+P338+P310: IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a doctor. P391: Collect spillage. **Storage.** P210: Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. P220: Keep away from clothing and other combustible materials. Store in a well-ventilated place. Keep container tightly closed. Keep dry.

Efficacy

Microorganism	Time
Bactericidal	**60 sec**
Pseudomonas aeruginosa	60 sec
Staphyloccocus aureus	60 sec
Enterococcus hirae	60 sec
Escherichia coli	60 sec
MRSA	60 sec
Yeasticidal	**60 sec**
Candida albicans	60 sec
Enveloped viruses	**60 sec**
Hepatitis B virus (HBV)	60 sec
Hepatitis C virus (HCV)	60 sec
HIV	60 sec
Measles virus	60 sec
Hepatitis Delta virus	60 sec
Influenza virus	60 sec
Coronavirus	60 sec
Herpesviridae	60 sec
Filoviridae (Ebola)	60 sec
Rubella Virus	60 sec
Rabies Virus	60 sec
Poxviridae	60 sec
Non enveloped viruses	**60 sec**
Adenovirus 5	60 sec
Murine norovirus	60 sec

UK	**Danger.**
BG	**Опасно.**
CS	**Nebezpečí.**
DA	**Fare.**
DE	**Gefahr.**
EL	**Κίνδυνος.**
ES	**Peligro.**
ET	**Ettevaatust.**
FI	**Vaara.**
FR	**Danger.**
HR	**Opasnost.**
HU	**Veszély.**
IT	**Pericolo.**
LT	**Pavojus.**
LV	**Bīstami.**
NL	**Gevaar.**
NO	**Fare.**
PL	**Niebezpieczeństwo.**
PT	**Perigo.**
RO	**Pericol.**
SK	**Nebezpečenstvo.**
SL	**Nevarno.**
SV	**Fara.**

CE 1639

UK CA 0120

Manufactured in the UK by:
Topdental (Products) Ltd.
Unit 1, Holmfield Industrial Estate,
Halifax, W. Yorks, HX2 9TN. UK
Tel +44 (0)1535 652750
www.topdental.co.uk

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