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
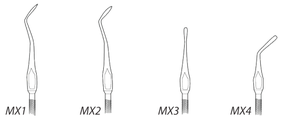
EN Reusable dental tools for handling and application of prepared measure of MTA material

BEFORE USE PLEASE READ THE FOLLOWING INSTRUCTIONS CAREFULLY

TO BE USED BY DENTIST ONLY

NON - STERILE PRODUCT. STERILIZE DIRECTLY BEFORE USE

The tools are made of corrosive – resistant stainless steel.

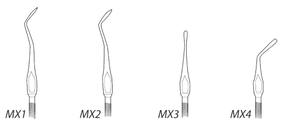HANDLING DENTAL TOOLS MADE OF STAINLESS STEEL
Special processes have been undertaken during the manufacture of the dental instruments to ensure corrosion-resistance. In spite of this, the durability of these dental instruments depends on their proper use and correct maintenance.

CLEANING, RINSING AND WASHING
The packaged instruments are non – sterile. Before sterilization, these new instruments should be washed thoroughly with warm water and using detergent. Before each use these dental instruments must be steriilzed.

STERILIZATION
All (reusable) dental instruments must be disinfected, then washed out to remove detergent and then dried before sterilization. Ultrasonic cleansing is recommended in particular as the most effective sterilization. Instruments should be sterilized in accordance with national requirements.

Recommended sterilisation methods:

Method	Minimal parameters:
Steam autoclave	134°C, for 5 minutes

PRECAUTIONS
Adhere to the guidelines for cleaning and sterilization of the dental instruments.
Any serious incident involving the product should be reported to the manufacturer and the competent authority of the State in which the user or patient is domiciled.

STORAGE
Dental instruments should be stored in temperatures ranging between 5-30°C and relative humidity below 70%.
Before use protect the instrument against contamination or damage of the material structure (cracks, bends, fractures, peeling, etc.). Failure to adhere to the above stated instructions decreases instrument durability and reduces
instrument lifetime.

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